

medical **products** research



# INSPIRATION IN RESPIRATORY THERAPY





## FREE FLOW

Active humidifier with a high flow generator  
and expiratory flow accelerator  
(EFA®).

# FREE FLOW


**FREE FLOW** allows the administration of high flow air, previously heated and humidified, if necessary oxygen-enriched, for an effective and comfortable ventilatory support, associated with a deep and efficacious airways clearance.

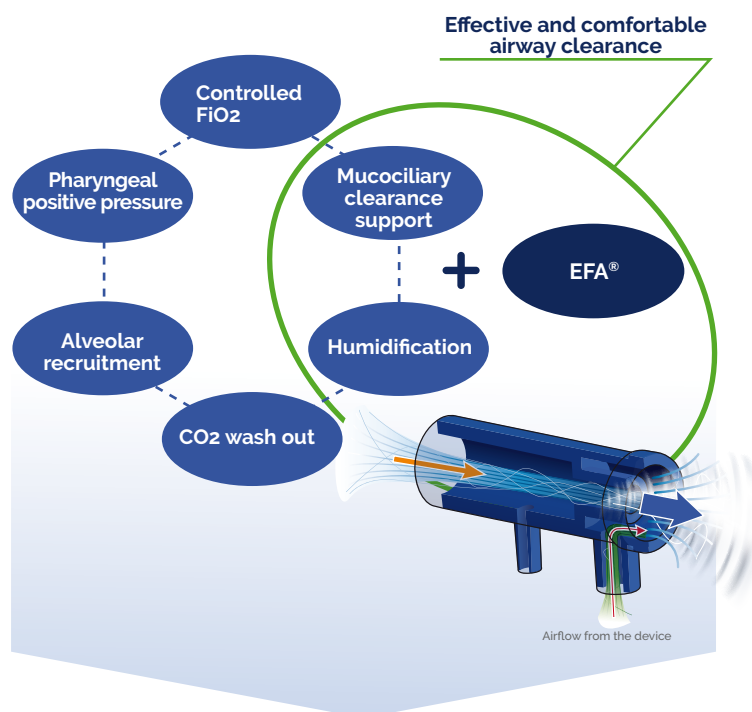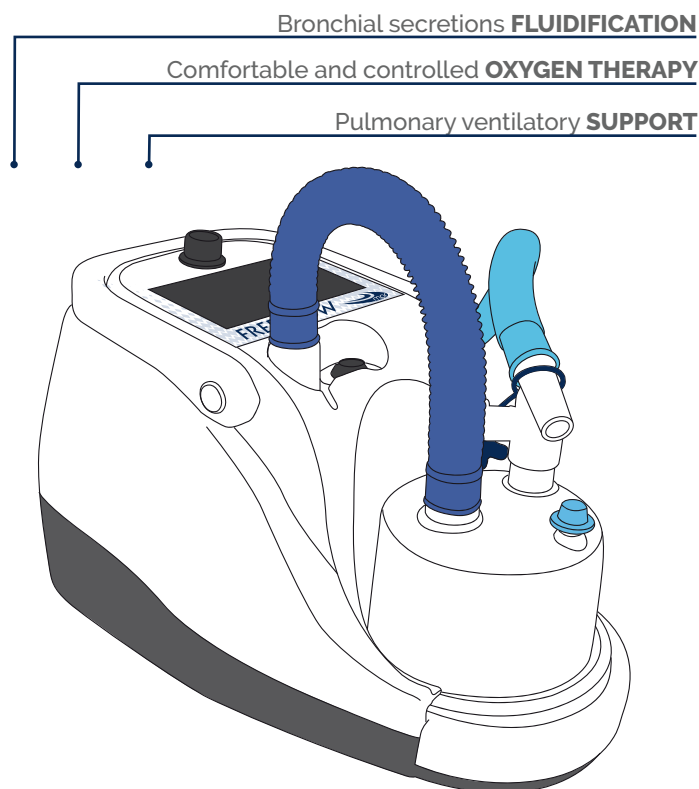

Besides the humidification therapy, the expiratory flow accelerator (EFA®) enables mucociliary clearance when the latter is compromised or ineffective, by facilitating the mobilisation and removal of secretions in a non-invasive way.

## FEATURES

- **FLOW SETTING** (2-60 L/min)
- **TEMPERATURE SETTING** (31-37°C)
- **CONSTANT FiO2 MONITORING**
- **INTEGRATED PULSE OXIMETER**  
Monitoring of SpO2 and heart rate
- **ELECTRONIC TEMPERATURE CONTROL**
- **SAFETY ALARMS**
- **5 LEVELS** of expiratory flow acceleration (EFA®)

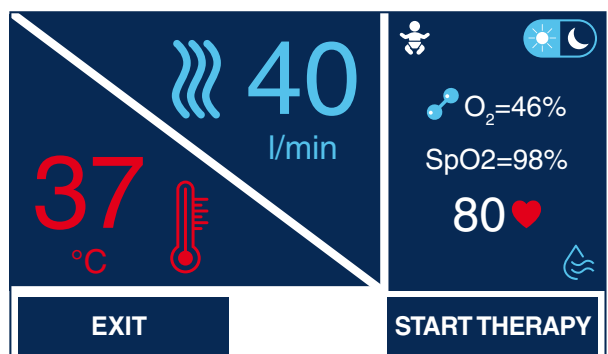

## VERSATILITY

- Vast range of **INTERFACES**
- **MODULABLE** accessories
- **DAY AND NIGHT** therapy
- **ADULT AND PAEDIATRIC** patient

## Technical Data

Power supply:	240V – 50Hz
Electrical absorption:	0.7 A
Electrical protection class:	II
Applied part type:	BF
IP protection grade (bag):	IPX1
Dimension and weight:	20.5 x 22 x 15 cm – 3 Kg

CE 1936

Medical Products Research Srl states that FREE FLOW is compliant with the DM 93/42/CEE directive, class IIa.  
The safety of the device is verified in accordance with the issued international standards.

This technical specification is updated at the printing time and may change according to our policy of continuous improvement. ©MPR srl

Information intended for healthcare professionals

**Medical Products Research S.r.l.**  
Via Novara, 68 - 20025 Legnano (MI)  
Tel +39 0331 597 992 - Fax +39 0331 485 089  
info@mpr-italy.it - www.mpr-italy.it

